

Evacuation Chair

for Public Access

Evacuation Chair

The comprehensive evacuation solution

Emergencies arrive unexpectedly. In that moment, preparedness is key. Are you prepared to evacuate occupants from your facility when facing an emergency?

Stryker's Evacuation Chair gives you the ability to evacuate disabled occupants weighing up to 500 pounds (with two operators), from multilevel facilities in emergency situations.

The patented Stair-TREAD[™] system assists in an easy descent down stairs by a single operator, with passengers weighing up to 200 pounds. Our extendable head and foot end lift handle design supports proper ergonomic lifting technique, while improving line of vision to feet and steps when lifting the chair on stairs or over obstacles.

Restraints

with Fastex[®] buckles for passenger security


Contoured seat for comfort

Features and benefits

- Assists in evacuation of disabled occupants from multilevel facilities
- Stair-TREAD system enables one person to control descent down stairs for passengers weighing 200 pounds or less
- Helps comply with fire safety requirements and emergency response guidelines
- Extendable head and foot end lift handle design supports proper ergonomic lifting technique, while improving line of vision to feet and steps when lifting the chair on the stairs or over obstacles
- Designed to maximize mobility in tight spaces
- Folds to compact size for convenient storage

Foot end lift handles extend for ergonomic lifting

The Evacuation Chair folds and stores in compact spaces.

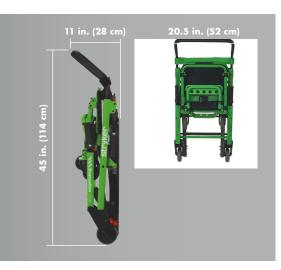

Foldaway foot support provides additional passenger security and stability

Upper control handle

extends to be ergonomically positioned for optimal operator control



Optional accessories





Wall mount 6254-050-000

Foot support 6254-029-000





In-service video 6253-001-167

Alarm 6254-045-000



Chair cover 6254-055-000

Storage cabinet 6254-002-000



Aggressive Stair-TREAD belt 6254-070-000



6250-140-000

Allows operator to control movement down stairs without lifting which may help to reduce the risk of lower back disorder.¹ The dual treads, made with Kevlar[®], stand up to tough use

Optional Aggressive Stair-TREAD

The Aggressive Stair-TREAD belt is designed to improve the grip of the Stair-TREAD system on steps that have a smooth surface with a rounded edge. To determine if the Aggressive Stair-TREAD belt is appropriate for your application, please contact your local Stryker sales representative.







Specifications

Evacuation Chair

| 6254-000-000 |
|---------------------|
| 45 inches (114 cm) |
| 20.5 inches (52 cm) |
| 28 inches (71 cm) |
| 11 inches (28 cm) |
| 34 pounds (15.5 kg) |
| 500 pounds (227 kg) |
| |

Foot support

| Foot support option part number | 6254-029-000 |
|------------------------------------|---------------------------|
| Step height from floor | 4.5 inches (11.4 cm) |
| Floor clearance | 2.25 inches (5.7 cm) |
| Width | 12.25 inches (31.1 cm) |
| Depth | 3.5 inches (8.9 cm) |

Evacuation chair cabinet

| Model number | 6254-002-000 |
|--------------|----------------------|
| Height | 47 inches (119.4 cm) |
| Width | 24 inches (61 cm) |
| Depth | 10.25 inches (26 cm) |
| Weight | 42 pounds (19 kg) |
| | |

Reference

1. Fredericks TK, Choi SD, Butt SE, et al. Postural Analysis of Paramedics Using Stairchairs. Human Performance Institute, Department of Industrial and Manufacturing Engineering Western Michigan University, Kalamazoo, MI 49008-5061, USA

Contact your local sales representative or visit stryker.com

Emergency Care Public Access

The information presented is intended to demonstrate Stryker's product offerings. Refer to operating instructions for complete directions for use indications, contraindications, warnings, cautions, and potential adverse events, before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area. Specifications subject to change without notice.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: Stair-TREAD, Stryker. All other trademarks are trademarks of their respective owners or holders.

The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.



Manufactured by:

3800 E. Centre Avenue Portage, MI 49002 U.S.A. P + 1 269 329 2100 Toll free + 1 800 327 0770 stryker.com

Distributed in Canada by:

Stryker Canada 2 Medicorum Place Waterdown, Ontario L8B 1W2 Canada Toll free 800 668 8323

Mkt Lit-2010 24 JUN 2020 REV B Copyright © 2022 Stryker